


Model 480 Fluxmeter RoHS



Model 480 features

- 5¾-digit DC resolution (1 part out of $\pm 300,000$)
- Automatic drift compensation
- Very fast peak capture
- AC frequency response to 50 kHz
- IEEE-488 and serial interfaces
- Storage of parameters for up to 10 existing coils
- CE mark certification

Model 480 Fluxmeter

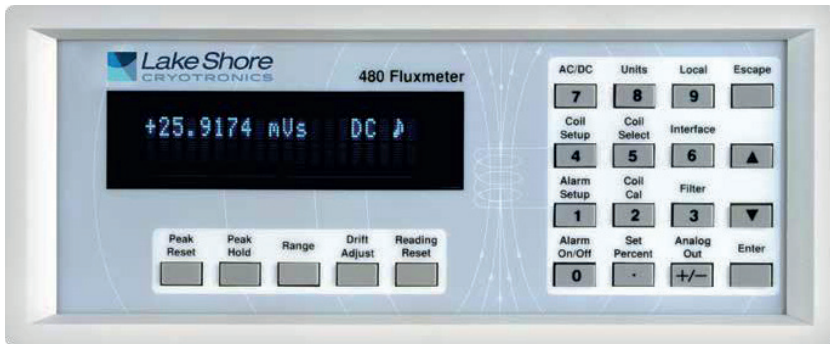

Product description

An advanced tool designed primarily for use in industrial and measurement systems settings, the Model 480 fluxmeter measures total flux from which B, flux density, and H, magnetic field strength, can be determined. The Model 480 is valuable for magnetizing, manual and automated magnet testing and sorting, and as the main component in BH loop or hysteresis measurement system applications. The Model 480 is compatible with most sensing coils and fixtures.

Manual magnet testing

A bright display and fast update rate make the Model 480 ideal for manual magnet sorting and testing. The instrument's low drift improves productivity with fewer adjustments. Remote terminals allow for foot pedal reading reset to keep hands on the work, not the instrument. Configurable alarms give an audible signal or relay closure to signify pass/fail.

Automated magnet testing

In automated testing, time is money. The Model 480 has many features to enhance throughput. The instrument has a fast update rate and fast settling time. It recovers quickly from reading reset to start a new reading cycle. The IEEE-488 and serial interfaces included with the Model 480 can be used to control most instrument functions. Relays and analog outputs can be used for automation without a computer interface.

Magnetizing

The magnetizing process places unique demands on all associated electronics. The

Model 480 responds with a very fast peak capture that can keep up with the fastest magnetizing pulses. Both a positive and negative peak can be captured from the same pulse. The input of the Model 480 is protected against the high voltages present during magnetizing.

Materials analysis

High resolution and low drift define a fluxmeter's role in analytical measurement. The high resolution of the Model 480 is reinforced by a low noise floor. A configurable filter helps keep the readings quiet. Automatic and manual drift adjustment modes help optimize the integrators' low drift characteristics. The IEEE-488 and serial computer interfaces included with the 480 allow automated data taking.

AC magnetic fields

Sensing coils are sensitive to AC magnetic fields but many conventional integrating fluxmeters can not measure AC fields. The Model 480 has an AC mode that enables it to measure fields over a wide frequency range using simple sensing coils. Applications are limited to field volumes as large as or larger than the coil, but for some it is an inexpensive way to make low drift AC field measurements.

Drift adjustment

Adjusting or nulling the drift of an analog integrator wastes time—it can be the only unpleasant part of using an integrating fluxmeter. Lake Shore innovation brings some relief. The Model 480 has a built-in drift algorithm that continually adjusts drift when

the instrument and coil are idle. It is ready when you are to make precision low-drift measurements. The adjustment algorithm has no effect during flux integration. Manual drift adjustment is also available.

Display

The Model 480 has a 2-line by 20-character vacuum fluorescent display. During normal operation, the display is used to report field readings and give results of other features such as max/min or relative. When setting instrument parameters, the display gives the operator meaningful prompts and feedback to simplify operation. The operator can also control display brightness.

Following are three examples of the various display configurations:



Normal reading—the default mode with the display of the live DC flux reading



DC positive and negative peak on—the display shows both the positive and negative DC peak readings



Alarm on—the alarm gives an audible and visual indication of when the flux value is selectively outside or inside a user-specified range

Model 480 Fluxmeter

FH-series Helmholtz coils

Lake Shore coils can be used with the Model 480 fluxmeter as well as with other fluxmeters. When used with a Model 480 fluxmeter, calibration and set up data are automatically loaded into the instrument. These probes and coils are accurately calibrated, using field standards maintained at Lake Shore. Most standards are traceable to physical standards such as a coil or probe of carefully controlled dimensions, or in some cases, to proton resonance. The coil constants are measured on the basis of the field generated by a current through the coil.




See pages 56, 57, and 58 for more information about available Helmholtz and search coils.

Helmholtz and search coils

Coils and probes wound by the user or from other manufacturers can be easily used with the Model 480. The Model 480 allows the user to save parameters for up to 10 existing coils/probes and quickly switch between them. Lake Shore also offers several sensing coils and probe assemblies for use with the Model 480 that have several conveniences. They are factory calibrated for accuracy and interchangeability. Calibration data is loaded into memory in the probe connector so it does not have to be entered by the user. Special coil assemblies are also available and can be designed to meet customer specifications.

Model 480 rear panel

- 1 Line input assembly
- 2 Serial I/O interface
- 3 IEEE-488 interface
- 4 Terminal block (for relays and analog signals)*
- 5 Coil input for user coils
- 6 Probe input for Lake Shore probes

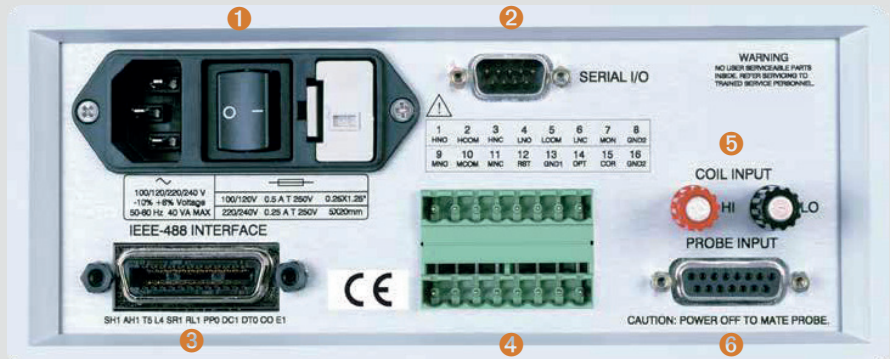

* The Model 480 terminal block has connections for external reset. With this feature, a foot pedal or programmable logic controller (PLC) can be used to start a new measurement cycle. The external reset is TTL-compatible and a logic low will activate a reset. The signal is internally pulled up to allow operation with a simple switch closure between pins 12 and 13.

Model 480 Fluxmeter

Model 480 specifications

Measurement

Number of inputs: 1

Input type: 2-lead, ground referenced

Input resistance: 100 kΩ or 10 kΩ

Maximum operating input voltage: 60 V

Absolute maximum input voltage: 100 V—**WARNING**—voltages between 60 V and 100 V will not damage the instrument but could result in personal injury or damage to other instruments

Update rate: 5 rdg/s display; 30 rdg/s IEEE-488; 30 rdg/s serial

DC

DC resolution: To 5¾ digits

DC integrator capacitance: 1 μF nominal

DC input resistance	100 kΩ	100 kΩ	10 kΩ	10 kΩ
DC ranges	300 mVs	30 mVs	30 mVs	3 mVs
DC resolution	0.001 mVs	0.0005 mVs	0.0005 mVs	0.0005 mVs

DC accuracy: Offset: ±10 μVs ±DC integrator drift; gain: ±0.25% of reading (<10 Vs/s maximum rate of change)

DC minimum dΦ/dt: 20 μVs/min

DC maximum dΦ/dt: 60 Vs/s

DC integrator drift: ±1 μVs/min, 0.0004% full scale/min on 300 mVs range (100 kΩ input resistance constant temperature environment)

DC peak

DC peak resolution: 4¾ digits

DC peak integrator capacitance: 1 μF nominal

DC peak input resistance	100 kΩ	100 kΩ	10 kΩ	10 kΩ
DC peak ranges	300 mVs	30 mVs	30 mVs	3 mVs
DC peak resolution	0.01 mVs	0.001 mVs	0.001 mVs	0.001 mVs
DC peak minimum reading	0.05 mVs	0.005 mVs	0.005 mVs	0.005 mVs

DC peak accuracy: Offset: ±100 μVs ±DC integrator drift; gain: ±5% of reading (<10 Vs/s maximum rate of change)

DC peak maximum dΦ/dt: 60 Vs/s

DC peak update rate: Reduces update rate to ¼ normal

AC

AC resolution: 4¾ digits (reduced to 3¾ digits on the 30 μVs range)

AC integrator capacitance: 0.1 μF nominal

AC input resistance	100 kΩ	100 kΩ	100 kΩ	100 kΩ
AC ranges	30 mVs	3 mVs	300 μVs	30 μVs
AC resolution	0.001 mVs	0.0001 mVs	0.01 μVs	0.01 μVs
AC minimum reading	3.000 mVs	0.3000 mVs	30.00 μVs	3.00 μVs

AC frequency response: 2 Hz to 50 kHz

AC accuracy: ±1% rdg ±10 μVs (10 Hz to 10 kHz sinusoidal); ±5% rdg ±10 μVs (2 Hz to 50 kHz sinusoidal)

AC integrator drift: NA

AC peak

AC peak resolution: 3¾ digits

AC peak integrator capacitance: 0.1 μF nominal

AC input resistance:	100 kΩ	100 kΩ	100 kΩ
AC peak ranges:	30 mVs	3 mVs	300 μVs
AC peak resolution:	0.01 mVs	0.001 mVs	1 μVs
AC minimum reading:	0.01 mVs	0.001 mVs	5 μVs

AC peak accuracy: ±5% rdg ±10 μVs (10 Hz to 10 kHz sinusoidal); ±10% rdg ±10 μVs (2 Hz to 50 kHz sinusoidal)

AC peak update rate: Reduces update rate to ¼ normal

Front panel

Display type: 2-line × 20-character vacuum fluorescent display

Display resolution: To ±5¾ digits

Display update rate: 5 rdg/s

Display units: Vs, MxN, WbN, VsΦ, MxΦ, WbΦ, G, T, Wbcm, A, %

Units multipliers: p, n, μ, m, k, M, G

Annunciators: AC—AC input signal, DC—DC input signal, Δ ∇ —positive and negative peaks, R—remote operation, \rightarrow —alarm on

Keypad: 21 full-travel keys

Interfaces

IEEE-488.2 capabilities: SH1, AH1, T5, L4, SR1, RL1, PPO, DC1, DTO, CO, E1

Serial interface: RS-232C electrical, DA-9 connector, 9600 baud

External reset type: Contact closure

Alarms

Number: 2

Settings: High and low setpoint, Inside/Outside, Audible

Actuators: Display annunciator, beeper, relays

Relays

Number: 3

Contacts: Normally open (NO), normally closed (NC), and common (C)

Contact rating: 30 VDC at 2 A

Operation: Follow high, low alarms with third relay indicating no alarm state—can be operated manually

Connector: Detachable terminal block

Monitor analog output

Scale: ±3 V = ±full scale on Vs range

Accuracy: ±1% of reading ±10 mV, (DC to 10 kHz); ±5% of reading ±10 mV, (10 kHz to 50 kHz)

Minimum load resistance: 1 kΩ

Connector: Detachable terminal block

Corrected analog output

Scale: User selected

Range: ±10 V

Resolution: 0.3 mV

Accuracy: ±2.5 mV

Minimum load resistance: 1 kΩ

Connector: Detachable terminal block

General

Ambient temperature: 15 °C to 35 °C at rated accuracy, 5 °C to 40 °C with reduced accuracy

Power requirement: 100, 120, 220, 240 VAC, +5% -10%, 50 or 60 Hz, 20 VA

Size: 216 mm W × 89 mm H × 318 mm D (8.5 in × 3.5 in × 12.5 in), half rack

Weight: 3 kg (6.6 lb)

Approval: CE mark, RoHS

Ordering information

Part number	Description
480	Model 480 fluxmeter

Please indicate your power/cord configuration:

1	100 V—U.S. cord (NEMA 5-15)
2	120 V—U.S. cord (NEMA 5-15)
3	220 V—Euro cord (CEE 717)
4	240 V—Euro cord (CEE 717)
5	240 V—U.K. cord (BS 1363)
6	240 V—Swiss cord (SEV 1011)
7	220 V—China cord (GB 1002)

Accessories included

106-739	Two 8-pin terminal block mating connectors
119-028	Model 480 user manual

Accessories available

4005	1 m (3.3 ft) long IEEE-488 (GPIB) computer interface cable assembly—includes extender required for simultaneous use of IEEE cable and auxiliary terminal block
CAL-480-CERT	Instrument recalibration with certificate
CAL-480-DATA	Instrument recalibration with certificate and data
CAL-NEW-DATA	Calibration data for a new instrument
RM-1/2	Rack mount kit for mounting one Model 480 in 483 mm (19 in) rack
RM-2	Rack mount kit for mounting two Model 480s in 483 mm (19 in) rack

Coils—see pages 56, 57, and 58 for more information

FNT-6R04-100	100 cm ² search coil
FNT-5P04-30	30 cm ² search coil
FH-2.5	Helmholtz coil, 64 mm (2.5 in) ID
FH-6	Helmholtz coil, 152 mm (6 in) ID
FH-12	Helmholtz coil, 305 mm (12 in) ID
FCBL-6	User programmable cable with PROM, 1.5 m (5 ft) long

All specifications are subject to change without notice