


Animal Monitoring and Gating Device

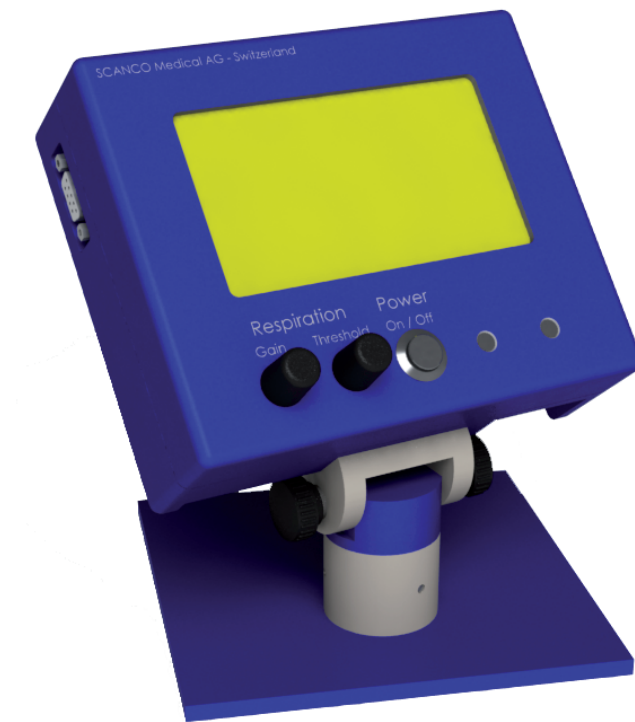

SCANCO MEDICAL

Animal Monitoring & Gating Device

The Scanco Medical Animal Monitoring and Gating Device allows users to

- Monitor respiration
- Monitor temperature
- Gate microCT scans (vivaCT 40 & 75) using the respiration signal

The unit includes

- Monitoring and gating device (main unit)
- Table stand
- vivaCT mount
- Respiration sensors
- Temperature sensor
- Adjustable gain and threshold for gating
- All connection cables/hoses



Specifications

Use recommended for	vivaCT 40 & 75
Main unit size	166 x 138 x 55 [mm], WxHxD
Main unit weight	1.3 kg
Power supply	100-240 V AC, 400 mA, IEC 60950
Respiration sensor	Pressure pad
Temperature sensor	PT100, resistance thermometer in stainless steel housing
Temperature range	5 ... 50° C (42 ... 122° F)



SCANCO Medical AG
Fabrikweg 2
8306 Brüttsellen
Switzerland
tel +41 44 805 98 00
fax +41 44 805 98 01

SCANCO USA, Inc.
985 Old Eagle School Road, Suite 511
Wayne, PA 19087
USA
tel 610 688 1440
fax 610 688 4976



Made in Switzerland