


Cooling/Heating Device

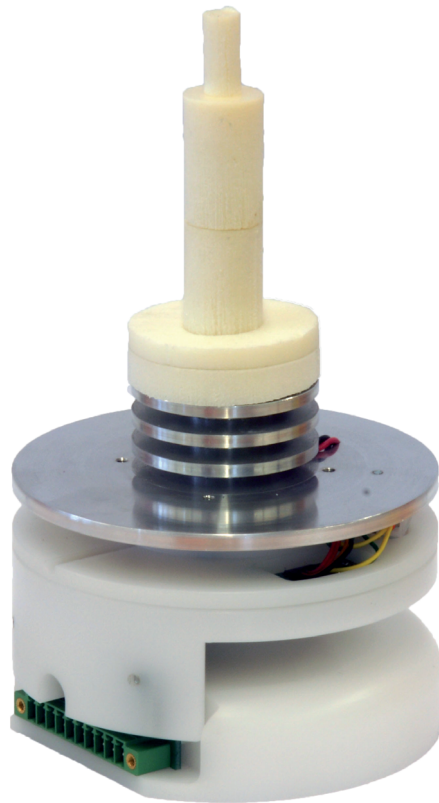


SCANCO MEDICAL

Cooling/Heating Stage

The SCANCO Medical Cooling/Heating Stage allows users to control the temperature down to 15°C below ambient temperature or up to 15°C above ambient temperature. It enables user controlled induction of temperature variations over time and monitoring and recording of set and actual temperature.

The unit includes

- Cooling/Heating stage
- Interface
- RS232 cable to connect interface with PC/Notebook
- USB to RS232 connector
- PC Software for temperature control

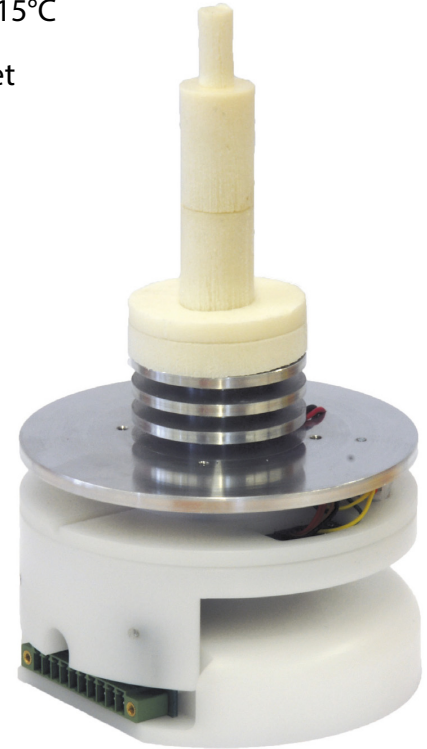

Specifications

- | | |
|-----------------------|---|
| - Use recommended for | μ CT 35, 40, 45, 50, 80, 90, 100, 100HE, vivaCT 40, 75 and 80 |
| - Max sample diameter | 8 mm |
| - Max sample height | 15 mm |
| - Accuracy | +/- 0.5°C |
| - Max cooling effect | -15°C from ambient temperature |
| - Max heating effect | +15°C from ambient temperature |



SCANCO Medical AG
Fabrikweg 2
8306 Brüttisellen
Switzerland
tel +41 44 805 98 00
fax +41 44 805 98 01

SCANCO USA, Inc.
P.O. Box 646
Southeastern, PA 19399
USA
tel 610 688 1440
fax 610 688 4976



Made in Switzerland